


pria™

BY

BLACK+DECKER™

HOME CARE COMPANION




QUICK START GUIDE

OKPRIA.COM

QUICK START GUIDE

WHAT'S IN THE BOX

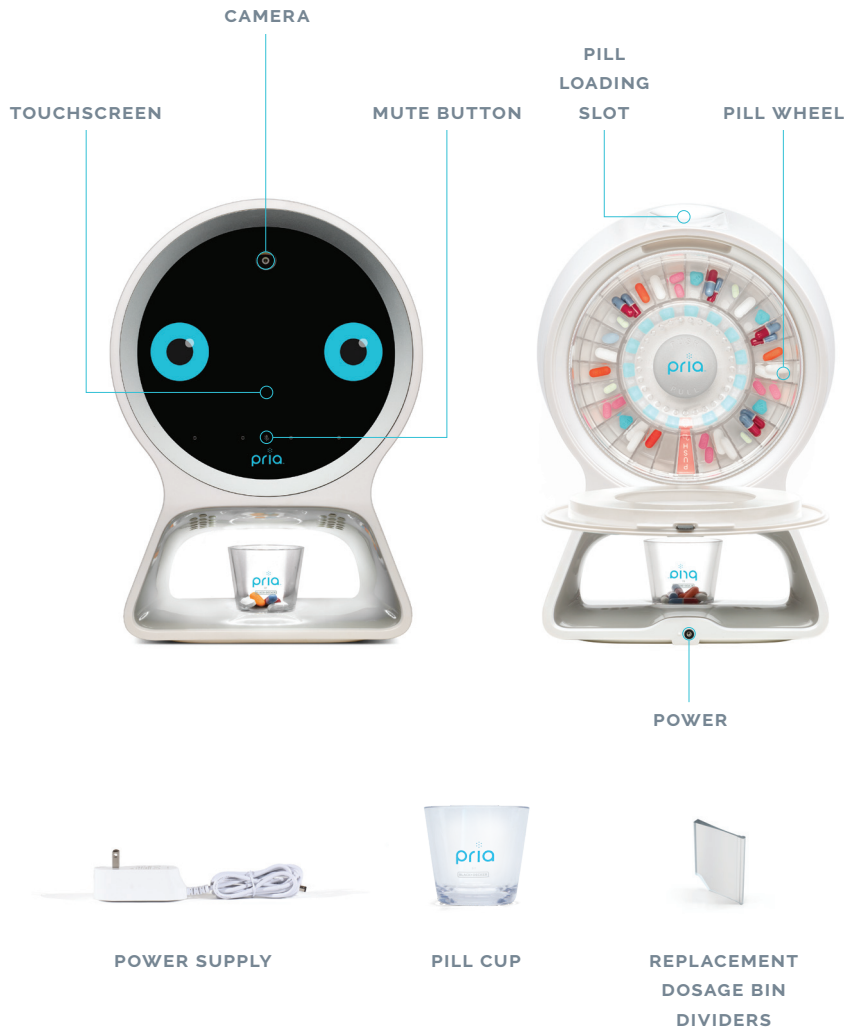

WHAT YOU WILL NEED



WiFi Network and Password



Your Medication



Your Smart Phone or Tablet
(not available for Windows devices)



Your Email Address

QUICK START GUIDE



FOLLOW THESE QUICK INSTRUCTIONS OR GO TO WWW.OKPRIA.COM/SUPPORT FOR A MORE IN-DEPTH LOOK.

STORAGE GUIDELINES

Place Pria on a flat surface where you usually keep your medications, in a dry area away from direct sunlight and heat.

IMPORTANT: Pria should not be stored in the bathroom.

1

POWER ON

Connect the power cable to the back of Pria and then to a power outlet. Pria's screen will light up and step-by-step instructions will appear, starting with setting up Pria's WiFi connection.



2

DOWNLOAD APP

Search for Pria By BLACK+DECKER in the Google Play™ or App Store® to install the app.

NOTE: The app should be installed on the phone or tablet of the person who will be using Pria. Your Pria contacts can install the app later on their own phone or tablet.



CONTINUED



3

PAIR APP WITH PRIA

Subscribe to your Pria plan to create your account. Tap "Go to settings" then "Profile" then select "Connect your Pria". Hold your device up and align it with the blue square on Pria's screen.

IMPORTANT: During onboarding, Pria will prompt you to enter a PIN. Write your PIN down and keep it in a safe place.



(Write your PIN here)

4

SETTING UP PRIA


Follow the on-screen instructions to take a few photos. Don't worry, these photos will only be used to safely dispense your medication only.

Listen to Pria describe her top features and benefits on her dashboard. Swipe from right to left with your finger to access the dashboard.

5

ADD YOUR MEDICATION SCHEDULE

To enter your medication schedule on the Pria app, tap the **ADD PILLS (+)** button on the top right corner. Follow the on-screen instructions and repeat for each medication.

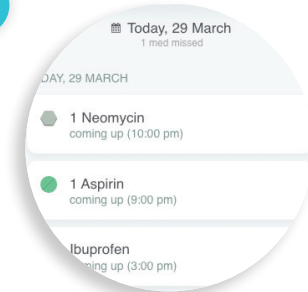

QUICK START GUIDE

6

LOADING MEDICATIONS

To start, say: **"OK PRIA. LOAD MY MEDICATIONS."** You can then choose between **ASSISTED** or **MANUAL FILL**.

ASSISTED FILL (PREFERRED METHOD)

With **ASSISTED FILL**, Pria will guide you with step-by-step instructions.

IMPORTANT: Be sure to add only one pill at a time as instructed. Don't worry about what dose is going in, Pria does that for you!

MANUAL FILL

If you opt for **MANUAL FILL** then follow the prompts to remove the Pill Wheel from the back of Pria. Once removed, follow the on-screen instructions.

To open each slot, push the red button with your thumb while rotating the clear protective cover with the other thumb.

IMPORTANT: When loading is complete, make sure the cover is in the **LOCK** position and place the wheel back inside Pria.

PLEASE NOTE:
The plastic dividers in the Pill Wheel will fall out if you turn the wheel upside down.



ALL DONE!

Any time you want Pria's attention, say "Ok Pria." When you see Pria's eyes appear on the screen, begin your request.



THINGS YOU CAN ASK PRIA

- "Ok Pria, volume five" (1-15)
- "Ok Pria, what's the weather this week?"
- "Ok Pria, how many calories are in a sweet potato?"
- "Ok Pria, what is my wellness plan?"
- "Ok Pria, how many tablespoons in a quarter cup?"
- "Ok Pria, dispense my meds"

"GOOD MORNING!
HOW DID YOU SLEEP LAST NIGHT?"



FOLLOW US:
@PriaAtHome





MODEL: PRIA-01A



INTENDED USE STATEMENT

The Pria By BLACK+DECKER Home Care Companion is intended for the medical purposes of:

a) providing alerts to the user, patients, or healthcare providers for pre-determined medication and nutritional supplement dosing schedules, and b) dispensing the appropriate oral dosage forms (pills, capsules, etc.) of solid medication and nutritional supplements.

SAFETY INFORMATION

1. Do not let power supply cord touch hot surfaces.
2. Do not mount near gas or electric heaters.
3. The use of accessory equipment not recommended by the manufacturer may cause an unsafe condition and void product warranty.
4. Do not use for other than intended purpose.
5. All servicing should be performed by qualified service personnel.
6. This device is not intended to be a diagnostic device.
7. Always check dispensed medication dose for accuracy.
8. Read all the information in the instruction manual and any other literature in the box before operating the unit.
9. To reduce the risk of injury keep hands off moving parts when in motion.
10. Wash Pill Wheel and Cup using water and soap, do not place in a dishwasher.

CAUTION

Do not overload power outlets. Plug the device into the appropriate voltage outlet.

WARNING

Do not remove battery.

Do not plug or unplug the power cord into the electrical outlet with wet hands.

Equipment should be mounted in locations and at heights where the camera has enough field of view to allow face detection and where the camera has enough light to operate.

FCC STATEMENT

This device complies with part 15 of the FCC rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. NOTE: The manufacturer is not responsible for any radio or TV interference caused by unauthorized modifications or changes to this equipment. Such modifications or changes could void the user's authority to operate the equipment. NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FCC CAUTION

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

IMPORTANT NOTE

When using the product, maintain a distance of 20 cm from the body to ensure compliance with RF exposure requirements.

INDUSTRY CANADA

This device complies with Industry Canada license-exempt RSS standard(s). Operation is Subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

Le présent appareil est conforme aux CNR d'Industrie Canada applicables aux appareils radio exempts de licence. L'exploitation est autorisée aux deux conditions suivantes: (1) l'appareil ne doit pas produire de brouillage, et (2) l'utilisateur de l'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement.

Pria By BLACK+DECKER
23 Drydock Ave, Suite 720W
Boston, MA 02210
Phone Number: 866-608-5168

FCC ID: 2AQK6S-EM010

IC ID: 24186-SEM010